

ReSound LiNX™

The world's smartest hearing aid

ReSound LiNX™ gives you smart new ways to fit – and satisfy – just about anyone.

The new frequency compression algorithm, Sound Shaper™, lets you fit ReSound LiNX on patients who need the flexibility of a fully featured BTE including telecoil and DAI functionality.

All of this in a 13 battery BTE that streams stereo sound directly from an iPhone®, iPad® or iPod touch® with ReSound LiNX 9 and 7 devices.

Options

- Available in 12 colors
- Direct Audio Input (standard feature but requires optional audio shoe)
- Thin Tube Adaptor for Open fittings
- Child Ear Hook
- Colored Side Indicator

Fitting Software

- Aventa fitting software version 3.7 or higher
- Wireless fitting with Airlink™
- Speedlink, HI-Pro, or NOAHlink programming interfaces




LN977-DW
LN777-DW
LN577-DW

ReSound LiNX 9
ReSound LiNX 7
ReSound LiNX 5

	ReSound LiNX 9	ReSound LiNX 7	ReSound LiNX 5
ReSound SmartRange chip	●	●	●
Surround Sound by ReSound			
Model			
WARP™ Compression - number of bands	17	17	9
Environmental Classifier	●	●	●
Sound Shaper™	●	●	●
Clean			
NoiseTracker™ II	●	○	○
Personalized Noise Reduction	●	○	○
WindGuard™	●	○	○
Expansion	●	○	○
Balance			
Binaural Directionality™	●	●	●
Directional Mix Processor	●	●	●
- Adjustable directional mix	●	●	●
Natural Directionality™ II	●	●	●
Synchronized SoftSwitching™	●	●	●
SoftSwitching™			●
AutoScope Adaptive Directionality™	●	●	●
MultiScope Adaptive Directionality™	●	●	●
Adaptive Directionality			●
Fixed Directionality	●	●	●
Binaural Environmental Optimizer™ II	●	●	●
Environmental Optimizer™		●	●
Stabilize			
DFS Ultra™ II	●	●	●
- Music Mode™	●	●	●
Auto DFS™	●	●	●
Convenience Features			
Synchronized Push Button	●	●	●
Synchronized Volume Control	●	●	●
SmartStart™	●	●	●
PhoneNow™	●	●	●
Comfort Phone™	●	●	●
iSolate™ nanotech	●	●	●
Fitting Features			
Gain Handles	9	7	6
Fully Flexible Programs	4	4	4
Synchronized Acceptance Manager	●	●	●
Onboard Analyzer™ II	●	●	●
In Situ Audiometry	●	●	●
Wireless connectivity			
2.4 GHz Wireless Technology	●	●	●
2.4 GHz Device-to-Device Communication	●	○	○
Direct audio streaming (Made for iPhone)	●	●	●
Wireless fitting with Airlink™	●	●	●
ReSound Unite™ TV Streamer 2	●	●	●
ReSound Unite™ Remote Control 2	●	●	●
ReSound Unite™ Phone Clip+	●	●	●
ReSound Unite™ Mini Microphone	●	●	●
ReSound Smart™ app	●	●	●
ReSound Control™ app	●	●	●

○ Standard
● Advanced
● Ultimate

Fitting Range

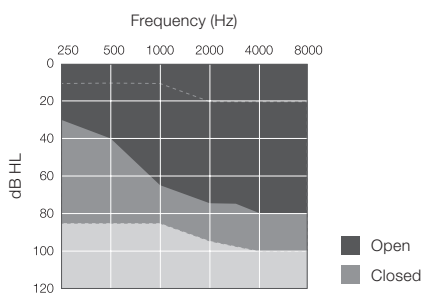

ReSound LiNX is compatible with iPhone 5s, iPhone 5c, iPhone 5, iPad Air, iPad (4th generation), iPad mini with Retina display, iPad mini and iPod touch (5th generation) using iOS 7.X or later. Apple, the Apple logo, iPhone, iPad and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

ReSound North America
8001 Bloomington Freeway
Bloomington, MN 55420
1-800-248-4327
customerexperience@gresound.com

ReSound Government Services
8001 Bloomington Freeway
Bloomington, MN 55420
1-800-392-9932
VAinfo@gresound.com

ReSound Canada
303 Supertest Road
Toronto, Ontario M3J 2M4
1-888-737-6863

ReSound

rediscover hearing

ReSound LiNX™


Technical specifications

		Open		Closed		
		IEC 60118-0 IEC 711 Ear simulator	IEC 60118-7 ANSI S3.22 2cc coupler	IEC 60118-0 IEC 711 Ear simulator	IEC 60118-7 ANSI S3.22 2cc coupler	
Reference test gain (60 dB SPL input)	1600 Hz/HFA	44	41	47	47	dB
Full-on gain (50 dB SPL input)	Max.	62	51	66	59	dB
	1600 Hz/HFA	53	46	56	52	
Maximum output (90 dB SPL input)	Max.	130	122	134	128	dB SPL
	1600 Hz/HFA	126	118	130	123	
Total harmonic distortion	500Hz	0.2	0.1	0.5	0.4	%
	800 Hz	0.5	0.2	0.6	0.4	
	1600 Hz	0.4	0.5	0.6	0.5	
Telecoil sensitivity (1 mA/m input)	Max.	91	-	94	-	dB SPL
Full-on Telecoil sensitivity @ 1mA/m	1600 Hz / HFA	81	75	84	81	
HFA – SPLIV @ 31.6 mA/m (ANSI)	HFA	-	99	-	106	
Equivalent input noise		23	23	23	22	dB SPL
Frequency range (DIN 45605/ANSI)		100-6900	100-6770	100-7080	100-6850	Hz
Current drain		1.1	1.2	1.1	1.2	mA


Data in accordance with IEC 60118-0, IEC 60118-7; Supply Voltage 1.3 V.

Patents pending


Maximum Output (OSPL 90)

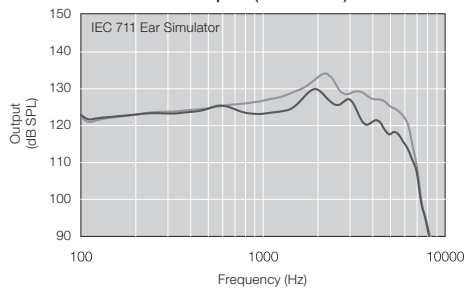


Maximum Output (OSPL 90)

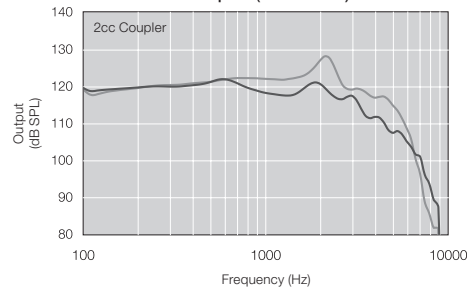


Full-On and Reference Test Gain

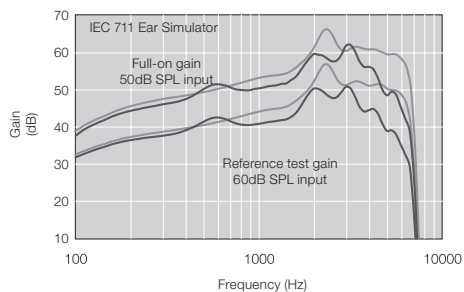


Full-On and Reference Test Gain

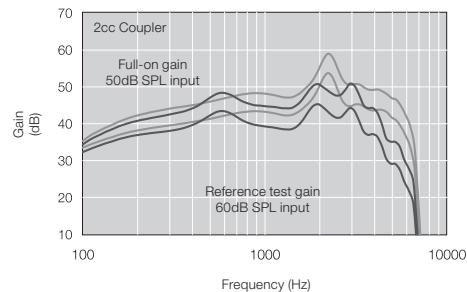

Telecoil response
Input level 10 mA/m

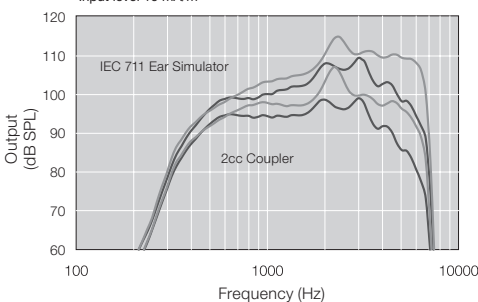

Input/Output Response

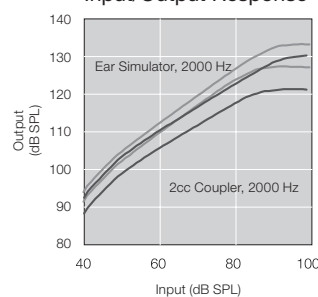

- Open
- Closed

Basic settings:

Full-on Gain, Reference Test Gain
MPO = Maximum Power Output
Maximum Band Width

Measured according to IEC 60 118-0 1983, amendment 1994; at 1.3 V, impedance 6.2 ohms and 23°C on O.E.S. according to IEC711 1981, resp on 2cc according to IEC60118-7 2nd edition 2005 (DIN average calculated at 500 Hz, 1000 Hz and 2000 Hz; HFA average calculated at 1000 Hz, 1600 Hz and 2500 Hz; 0 dB SPL sound pressure equals 20µPa). All measurements without DSP features activated unless indicated otherwise.

All specifications are subject to change without notice

400304011-US:14.07-Rev.A